

Extraordinary energy meter

Product instruction manual



MXSELLA

The uniqueness of the extraordinary energy meter

Extraordinary energy meter by Korea Mixiu Medical Co., LTD independent research and design. "RFPEMF is the fifth element of health" in the book, in youth, the human body emits pulse energy about 3-5 times per minute, this pulse energy can improve blood circulation, is the function of all parts can be absorbed fresh blood. RFPEMF- Pulsed electromagnetic energy therapy (PF Pulsed Electric Magnetic Field), belongs to the field of physical medicine treatment, PEMF is to simulate the earth's natural pulse magnetic electricity to improve the body's circulation metabolism. RFPEMF extraordinary energy meter is based on the German microcrystalline magnetic technology, citing Israel sine wave magnetic energy technology, the use of international advanced circuits, humanized design, exclusive patented technology, solemnly launched a can enhance the human body independent life vitality of the health instrument.

The principle of extraordinary energy meter

Extraordinary energy meter is sublimated on the basis of German microcrystalline magnetic technology, the use of international advanced circuits, humanized design, exclusive patented technology, through the sole microcrystalline magnetic resonance of the sinusoidal wave magnetic energy to the whole body, in the human cell polarization heat phenomenon, the positive and negative ions in the human body, the molecules that constitute the organization will produce violent vibration, rotation, reversal, collision and other movements, More than 60 points from the six meridians on the sole of the foot to the deep introduction of major meridians, dissolve blood stasis, dredge qi and blood, accelerate the body's microcirculation system, make the body instant fever, improve their own self-healing ability and independent life vitality of the rehabilitation health care instrument, is currently a hot spot in the health care and medical field. Thermal phenomena, positive and negative ions in the human body, constitute tissue.

Extraordinary energy meter technology

After the extraordinary energy is introduced into the human body, the body temperature rises rapidly, causing cell resonance, forming the thermal circulation and molecular polarization effect of the whole body, the molecular structure of blood and body fluids becomes smaller, the flow of qi and blood is accelerated, and the mucous membrane, old waste, cholesterol, fat and all harmful substances attached to the body fluids, blood and blood vessel walls are gradually discharged. The blood flow of capillaries is 4-5 times the usual, the vascular space is increased, the blood flow is increased, the blood carries oxygen and nutrients are greatly increased and reach the body smoothly, antibodies, cells and genes are repaired and restored to vitality; Activate the detoxification organs in the state of paralysis, and enhance the metabolic self-purification function and immunity of the self-repair system.



Applicable population

- 1.Obesity: Passively improve aerobic circulation and burn body fat.
- 2.Cardiovascular diseases:hypertension,hyperglycemia, hyperlipidemia, gout, varicose veins, etc.
- 3.sub-health symptoms: body deficiency qi cold, back pain, joint pain and other poor blood circulation.
- 4.Urological diseases: prostatitis, prostatic hyperplasia, etc.
- 5.gynecological diseases: palace cold, dysmenorrhea, irregular menstruation, adnexitis, cervicitis and so on.
- 6.immune diseases: cold, rhinitis, skin diseases, skin aging dry and cracked.



catalogue

1. Matters needing attention.....	1
1.1 Use contraindications and risk tips	1
1.2 Safety precautions	2
1.3 Precautions before use	3
1.4 Precautions for use	3
2. Introduction to use	4
2.1 Instrument introduction	4
2.2 Instructions for use	5
3. maintenance	6
3.1 Cleaning and maintenance	6
3.2 Q&A	6
4. Maintenance clause	8
5. Technical parameter	8

1. Matters needing attention

1.1 Use contraindications and risk tips

Because this instrument generates internal heat source and electromagnetic field, promotes blood circulation, and has a certain expansion effect on blood vessels, the following personnel are strictly prohibited to use:

- 1.Excessive drinking is strictly prohibited.
2. For patients with symptoms of high fever and bleeding, it is strictly prohibited to use it for those with wounds that have not healed after surgery.
3. Female menstrual period, infants, epilepsy and mental patients are strictly prohibited.
4. Do not use for children or people without the ability to express consciousness.
5. It is strictly prohibited to use during pregnancy, within 2 months of lactation and people with serious cardiopulmonary insufficiency.
6. Myocardial infarction, cerebral infarction patients, people suffering from cerebral haemorrhage and cerebral stroke, recovered less than half a year, the use of the instrument is prone to recurrence and risk (prone to cerebral haemorrhage and heart disease recurrence), so patients with the above diseases are strictly prohibited.
7. Electronic devices such as pacemakers, stents and cochlear implants, and metal objects such as steel plates and bone and joint stents implanted in the body are strictly prohibited.
8. Diabetic foot (ulceration, injury) is strictly prohibited. High blood sugar, numbness of the legs and feet, and insensitivity to heat should be used under the supervision of a professional, and the temperature should be lower than that of normal people (energy intensity should be adjusted from low to not too high) to prevent local burns caused by poor plantar contact.
9. People with weak constitution and heart and kidney failure are prone to dizziness and physical exhaustion, so it is not recommended to use.
10. Dizziness will occur when blood sugar is too low, so it is strictly forbidden to use this instrument on an empty stomach.
11. Minors (over 13 years old and under 18 years old) should be used as appropriate, and new users who do not know the instrument need to be supervised when using it.

1.2 Safety precautions

1. In the process of conditioning, if there is an emergency (such as power failure, failure, etc.), the power switch should be turned off in time.

2. Before operation, place both feet on the energy plate of the instrument, close to the energy plate, and then start the instrument. Pause the instrument before leaving the energy plate of the instrument with both feet. Note that during the conditioning process, the feet should not touch, so as not to affect the conditioning effect.

3. If you feel any abnormality of the device or body, please stop using it immediately.

4. Cut off the power supply in time after each use. Do not use the instrument for heating or cover the instrument with objects.

5. Do not use it in a place with high humidity such as the bathroom, or use it while taking a bath.

6. It can not be used with other electrotherapy instruments at the same time, which may cause burns and unstable output of the instrument.

7. Do not touch the energy plate area with metal objects.

8. When the equipment is in the working state, whether it is a maintenance test machine or a normal test machine, it is strictly prohibited to add power in the no-load state, otherwise the RF source is fully reflected, which will burn the high-frequency devices!

9. If smoke or odor occurs when the instrument is used, unplug the power supply and consult after-sales service.

1.3 Precautions before use


1. Before using this instrument, please remove the metal items on your body, such as rings, necklaces, watches and other items.
2. Drink more warm water when using the instrument (ensure that you drink 200 ml to 300 ml of warm water before, during and after the machine), or drink targeted health drinks.
3. Fasting, less than half an hour after a meal, half an hour after exercise and excessive fatigue should not be used immediately.
4. This instrument is an auxiliary product for sub-health people and sub-health rehabilitation and nursing, and shall not replace drugs and treatment. People suffering from diseases should be approved and guided by doctors.

1.4 Precautions for use

1. The feet should be completely close to the operating energy board when using.
2. Keep your feet and knees away from each other.
3. Do not touch the energy board with your fingers during use.
4. Do not touch your ankles with your hands to avoid static electricity.
5. If you feel uncomfortable during use, please lower the energy or stop using immediately.
6. During the conditioning process or within a few hours after the conditioning, the ankle is a normal reaction.
7. Users should not blow directly, and should avoid the machine at the vent.
8. After the machine, sweat should not be directly blown cold air, and it is not recommended to take a bath within 1 hour after use.
9. Adjust the energy intensity of the instrument according to your physical condition and acceptance during use, and carefully adjust the energy intensity according to the principle from low to high

2. Introduction to use

2.1 Instrument introduction

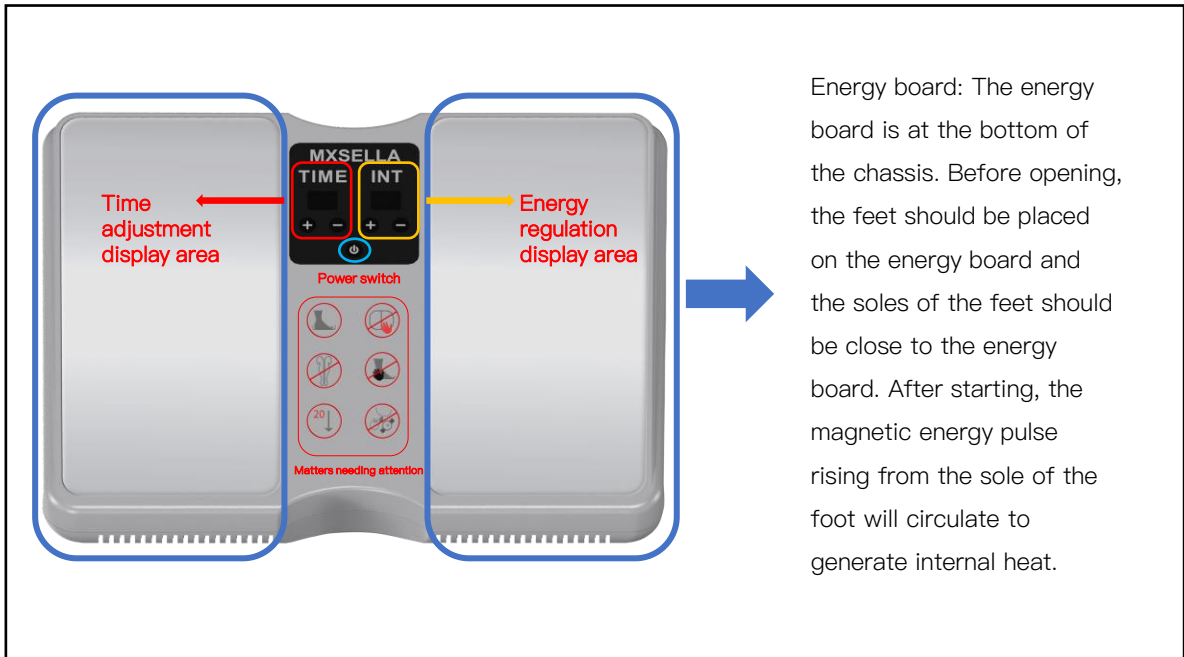

Use both feet fully close to the operating energy board



Do not touch your knees when using



If you feel uncomfortable during use, please lower the energy or stop using immediately



Do not touch the energy board with your hands during use



Do not touch your ankles with your hands to avoid static electricity



Please remove metal items such as rings, necklaces, watches and other items

2.2 Instructions for use

1. Before the machine, check whether the electrode plate cable is intact; It is recommended to remove metal items such as necklaces, rings and watches; It can be operated under the guidance of a person who has read the instructions carefully or knows how to use it.
2. In the ready state, connect the power cable, turn on the power switch, disinfect the energy board with 75% alcohol wipes, and place your feet on the energy board after the alcohol evaporates. Do not touch your knees, or touch the energy board and legs with your hands.
3. Use the instrument to start the function through the button, open the magnetic therapy, adjust the intensity of the magnetic energy pulse, you can start the leg and foot physiotherapy, during the process, the user according to their own tolerance to adjust the energy intensity, it is recommended to gradually from low to high, feel warm and comfortable. (Adjust the intensity of the energy output. Because the first user meridian stasis, energy concentrated in the ankle, easy to produce ankle pain, instep blood vessels tingling swelling feeling. The intensity of the energy of the instrument is adjusted according to the user's tolerance. At the beginning, the energy can be reduced. After 3 to 5 days, the time and energy can be adjusted according to the user's adaptation. When conditioning, drink more warm water to accelerate the circulation of qi and blood, and the discharge of toxins)
4. Turn off the instrument after the completion of use, the instrument in the stopped or suspended state (the working time is restored to 0 or stop the countdown), turn off the power switch, unplug the power cord, wipe the energy board with 75% alcohol, and store the instrument in a cool and dry place after the alcohol is volatilized.

3. maintenance

This equipment is composed of electronic instruments, and its maintenance and repair shall comply with the general regulations of electronic instruments. In storage and use, it is necessary to prevent bumping, extrusion and vibration.

3.1 Cleaning and maintenance

Before and after use, use 75% alcohol wipes to wipe and disinfect. If there are traces, wipe with a wet rag, then wipe and disinfect with 75% alcohol wipes, and then use after the alcohol evaporates. Store the instrument in a cool, dry place when not in use.

3.2 Q&A

1.ankle pain, instep blood vessels distension?

The ankle is where the meridians gather, especially near the three Yin Jiao points. When the Yang produced by the heat therapy instrument is following the meridians and collaterals, it is due to the stasis of the meridiansBlock, there will be ankle pain, muscle swelling feeling, when the Yang opens up this part, this sour swelling pain feeling will be reduced or even disappear.

2.why to drink warm water in the process of conditioning?

- A) Thermotherapy conditioning will make the body temperature rise, then the human body needs to replenish water, aerobic metabolism needs water and body fluids to participate in accelerated circulation.
- B) Deep heat therapy, the body surface sweat needs to replenish water, help the body to gather toxins with sweat discharge, so that detoxification more pleasant.
- C) The thermal effect of warm water helps to enhance the efficacy of hyperthermia.
- D) Add warm water to avoid the symptoms of hypoglycemia and dizziness in the weak.

3. How to deal with flustered discomfort during use?

If this happens, the use should be stopped immediately, the user should be provided with sugar and warm water, and allowed to rest properly, blood pressure measurement before leaving.

4. extraordinary energy meter how many times a day is appropriate?

It is recommended to use once a day, with treatment can be used twice, but need to be more than 3 hours apart.

5. What time of day is best to use?

It is best between 6:00 in the morning and 19:00 in the evening, if it is used too late, it is easy to cause Yang not to enter Yin, affecting the quality of sleep.

6. how many times to use for an adjustment cycle?

The early stage of health care can be 1 times a day, 10 days for a cycle, and then it can be adjusted once every 1–2 days. With chronic disease treatment and conditioning can be 30 days for a cycle, and then depending on whether the condition is adjusted to increase the conditioning cycle, such as basic improvement can be adjusted to adjust once every 1–2 days.

7. Why do I have sore knees or sore calves?

Knee soreness is the manifestation of wind cold and humidity emanating from the joint, and sore calf is the manifestation of muscle repair after simulating movement under high–frequency vibration.

8. why do you have lethargy or difficulty falling asleep, irritability?

Most of it is a phenomenon of physical deficiency, insufficient Yin blood, Yang energy, and short imbalance of nutrition that can not keep up with nutrition, and it is necessary to supplement nutrition and take appropriate health products or drugs.

4. Warranty clause

1.the main engine warranty for one year, accessories warranty for three months, lifelong maintenance!

- a) the main machine warranty for one yearThe main machine is guaranteed for one year, and the new machine can be replaced if there is non–human damage within one month.
- b) accessories warranty for three monthsParts warranty for three months, three months after the cost of parts repair fees.
- c) In case of man–made damage, the maintenance costs shall be borne by ourselves.

2.for the following personal reasons caused by the failure, do not enjoy the warranty service.

- a) Failure caused by unauthorized disassembly and modification of this product;
- b) Failure caused by careless fall and fall;
- c) failure caused by lack of reasonable maintenance;
- d) Failure caused by failure to follow the correct instructions in the instructions for use.

3.This product design specifications are subject to change without prior notice, if you have any questions, please contact the company.

5. Technical parameter

item	argument
technology	Biomimetic sine wave magnetic energy technology, constant temperature heating technology
Controls	Digital display, button touch
Operating voltage	220V
Output power	Static power: 20W; Maximum power: 350W

MXSELLA

Korea Mixiu Medical Co., LTD

www.mxsella.co.kr

Company address: No. 34, No. 342,
Suji Road, Suji District, Yongin, Gyeonggi Province